

Aircast® AirSelect™ Elite patient information

(01EP-XS, 01EP-S, 01EP-M, 01EP-L, 01EP-XL)

Aircast® AirSelect™ Elite

AirSelect Elite :

- Engineered to provide the ultimate in protection, comfort and oedema control.
- Specifically designed for stable fractures of the lower leg, foot and ankle; severe ankle sprains; and post-operative use.
- Pre-inflated front panel aircell and overlapping distal duplex aircells
- Enclosed toe cover.

Application Information


A cotton or cast sock will enhance comfort - INCLUDED

1 Prepare Boot


- Check aircells are deflated (see 'Deflate Aircells' below)
- Unfasten straps.
- Remove front panel.
- Open liner.

2 Apply Boot

- Put on the sock included with the boot. While seated, place foot in boot.
- Position heel against back of the boot.
- Wrap foot flaps over the foot and the toe cover around toes.

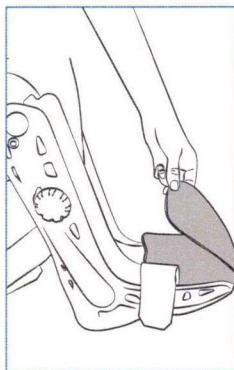



- Replace front panel.
- Position front panel so the Aircast logo is over the toes.
- Secure straps from bottom to top.
- Tighten until snug and comfortable.

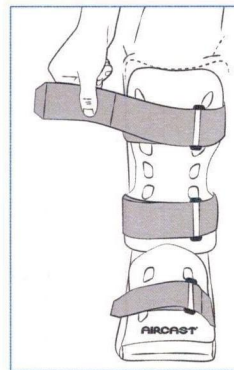

3 Adjust Aircell Compression


The ELITE walking brace aircells are initially deflated

- Begin inflation on injured side first. Turn the SELECTOR found on the right-side of the brace to the appropriate number. Number 1 inflates the aircell on the right-side, Number 2 inflates the left-side and Number 3 inflates the back aircell.
- To inflate aircell, press and release LARGE inflation pump found on the top right-side of brace, ensuring that the small hole is covered whilst doing so. Pump until aircell feels snug.
- Repeat to inflate other aircells, inflating Number 3 last.
- When all aircells are inflated, turn selector dial to "lock" position. This prevents accidental deflation.

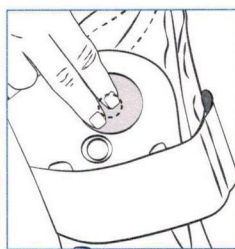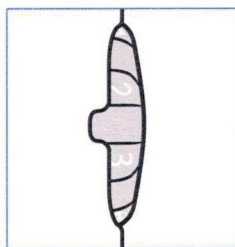

4 Deflate Aircells

- To deflate aircells, turn the SELECTOR to the appropriate aircell number.
- Press and hold finger on SMALL release button found on the top right-side of brace.
- Deflate until pressure is comfortable.

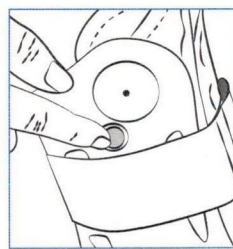

NOTE:

Deflate aircells before removing boot for proper fit when reapplying.

Higher pressure provides more support when walking, lower pressure is more comfortable when sitting or reclining. Pressure can also be adjusted by tightening or loosening straps.

If on an airline flight, deflate cells and loosen straps.

CLEANING INSTRUCTIONS:

The shell of the boot can be wiped with a damp cloth and a mild soap. Socks can be hand or machine washed up to 60 °C and tumble dried if required.

MATERIAL COMPONENTS

Shell: Polypropylene
Front Panel: High Density Polyethylene
Straps: Nylon
Lined Aircells (Anterior And Posterior): Polyurethane W/ Polyester Fabric
Inflation Pumps: TPU (Thermoplastic Polyurethane)
Liner: Nylon/Polyurethane/Polyester (Outer/Foam /Inner Mesh)
Binding: Polyester
Outer Sole / Insole: EVA/Rubber
Selector Dial: Nylon
Sock (Elite): Acrylic/Cotton

DJO Global
1a Guildford Business Park
Guildford
Surrey
GU2 8XG
UNITED KINGDOM
Tel: 0800 587 0857
www.DJOglobal.co.uk



DJOglobal.co.uk

© 2015 DJO, LLC 00-1319-EN Rev C